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CE

MOBIAK MEDICAL • GASES • FIREFIGHTING est. 1977

1. Introduction

We would like to thank you for choosing a MOBIAK SA stretcher. The product user manual is an important document that provides you with information on how to use it and technical features. Please read the user manual carefully before using the product. Learn the product and practice before any application. This manual applies to our model: **0805906.**

- The annotations and illustrations in this manual might be slightly different with the real parts due to quality improvement or changing design. Please check it carefully.
- The product you have in your hands may differ from the one in the following photos. The assembly process is similar.

SYMBOLS: Please refer on ANNEX A

2. Intended Use:

The product is used for the transportation of patients in emergency or in situations demanding a quick and safe evacuation of the patient. It has to be used exclusive by professional and qualified personnel, duly trained for emergency care. Can be used in ambulances and in hospital and sports environments.

3. CAUTION (!):

- In case of damage or malfunction of your product DOT NOT use the product and please contact the authorized and trained representative.
- Do not use the product for a purpose not indicated in this manual.
- MOBIAK S.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration on the frame of the product.
- MOBIAK S.A reserves the right to change the information contained in this document without previous notice.
- There is a ± 3%. Deviation for the Dimensions of the Products.
- It is always recommended to follow the procedures designed for emergency response.
- Improper use of the device can cause injury to the patient.

4. GENERAL WARNINGS:

- If you don't read this user manual, it is preferable not to use this device or another available part.
- If you don't understand the usage or the precautions, please contact with the dealer or the proper technical person before you use the device because may be caused damage.
- Carefully READ the entire manual before using the product.
- Save this document for future reference.
- The product should be used only by qualified personnel.
- Keep away the product from heat sources.
- Be careful when children are nearby and do not allow the children to play or operate with the device.
- Special care is always required where there are moving parts that can cause limb entrapment or injury.

5. Product Presentation:

- The stretcher is made of high-density polyethylene - HDPE and has aluminum tubes around the perimeter. This makes it robust, flexible and reliable for use in special situations.
- Its material is durable with high strength and non-absorbent. It is also easy to clean and disinfect.
- The maximum Load is 250Kg.
- It includes three (4) safety belts, with a car-style lock, 5cm width for tying up the injured and a lifting system with four lashing straps with safety rings, which fit into the special round holes on the stretcher.
- It has an adjustable footrest mechanism, which helps with both horizontal and vertical transport.
- Contains internally, an integrated soft, non-absorbent coating 1cm thick for the comfort of the wearer, which is easily cleaned and disinfected.
- It is foldable, you can divide it into two parts so that it can be stored and transported easily.
- It is suitable and functional for rescue needs either in the mountains or at sea.
- It is not compatible for use with X-ray and MRI devices. (X-Ray & MRI non-compatible)

6. Package Contents

The package contains: the main product, the adjustable foot fixing mechanism x1, the safety belt x 4, the lifting system x1 & the instructions for use x1.

7. Specifications

Reference Number	0805906
Name	Basket Stretcher
Frame	Foldable in two parts, for easy transport in vehicles
Frame Material	High density polyethylene – HDPE and aluminum tubes around the perimeter
Colour	Orange
Product Dimensions	218x62x16 cm
Product Dimensions Folded	124x62x20 cm
Maximum Load	250 Kg
Product Net Weight	20Kg with all accessories on and 15Kg without.
Safety Belts	4 pieces, car type, quick release, width 5cm
Lifting system	Yes, it is included in the package.
Adjustable foot fixing mechanism	Yes, it is included in the package.
Handles	12 peripheral handles & 4 round holes for placing safety rings for the lifting system.
Other	Internally it contains a soft, non-absorbent coating 1cm thick.
X-RAY & MRI Compatibility	NOT COMPATIBLE
Way of Transfer	Horizontal and Vertical
Master Carton Dimension / Gross Weight	130x65x23 cm / 1 pc per Master Carton / 23Kg

8. Before Each Use:

1. The actual need for use of this stretcher must therefore be carefully assessed by the rescuer in accordance with his/her training guidelines
2. Check if the frame is damaged in order to guarantee a safe use of the product (there are no cracks or fractures in the frame).
3. Check that the surface of the stretcher is not cracked.
4. In case of damage do not use the product and contact your dealer for further instructions.
5. Check, if the product is properly assembled and all screws are secure and well screwed on.
6. Always check the condition of wear for all the parts to ensure the absolute safe use of the product for people and objects.
7. Periodically check its sturdiness and the stability.

9. Use of the Product:

1. Join the two parts of the stretcher and secure with the pins.
2. Successfully move the patient or user onto the stretcher and secure with the belts. One belt should be in the chest area, one belt should be in the pelvic area and the other belt should be in the leg area to hold and securer these areas of the transferee.
3. Adjust the footrest, so it is in contact with the legs of the user.
4. The presence of 12 handles around the perimeter allows you to lift the patient by grasping the stretcher from various sides.
5. Before lifting the stretcher, make sure that the product in its folding parts is secured, you have properly strapped the patient with the safety belts, and the footrest fits the patient's legs.
6. DO NOT use with patients who exceed the weight limit on the product label and in this user manual. In specific cases or in the event of excessive patient weight, intervention by additional personnel may be required to ensure patient comfort and safety.
7. Always inform the patient when adjusting to the stretcher. Never leave stretcher unattended with patient inside

10. Requirements for operators:

- To be trained to deal with corresponding situations.
- Or have practical knowledge of handling patients in emergency situations.
- Be familiar with the procedures and the stretcher.
- Be physically able to assist a patient.

11. Maintenance:

Check the stretcher after each use. The inspection is intended to investigate anything that may compromise the integrity and proper operation of the stretcher. Check the product for any looseness, cracks, breaks, bent or damaged parts.

WARNING
In the event of a problem occurs, DO NOT use the product until the issue will be repaired.

12. Cleaning And Disinfection:

- Clean the plastic parts with water, rinse and wipe well with a soft and clean cloth and dry away from direct sunlight.
- You can use a wet cloth and mild soap on the surfaces that need a better cleaning.
- DO NOT leave the product wet or soiled with dust or mud.
- Do not use chemical cleaners to clean the frame. This may damage the product's coating surface.
- CAUTION: Do not use acidic or caustic solutions, **D**o not use solvents or stain removers
- If you must disinfect the product, use a common, mild disinfectant, after reading the instructions provided by the disinfectant manufacturer.
- Protect the product from scratches, cuts and punctures.

13. Conditions of Storage and Transport

The product should be stored packaged under the following conditions:

- For long-term storage the product should be placed in a dry and ventilated place and should not be placed in an environment with high and rapid changes of temperature.

- DO NOT store the product in locations that are close to high temperature sources or direct sunlight for long periods of time.
- The product should be isolated from chemicals and corrosive substances.
- The product should be packaged to protect it from dust, and parts that may be affected or damaged by oxidation should be lightly lubricated.
- During transportation and storage, when packed, the maximum stack height is six layers.
- The product can be transported by any common vehicle.
- Temperature range: -20°C ~ 60°C.

14. Environmental Protection Statement:

If one day you find that your product needs to be replaced or it is no longer working for you, consider protecting the environment:
1) Do not dispose your product along with the rest of the public waste (this is also the meaning of the shown recycling sign).
2) Contact your Public Authorities and they will instruct you of the Recycling centers to which your product must be disposed.
3) Correct disposal of your product helps the protection of the environment as well as the recycling of the product's components.

15. Incident Reporting:

Any serious incident that occurs in relation to use of the device must be reported to the manufacturer and a competent authority of the Member State in which the user and / or patient resides.

16. Declaration Of Conformity:

We are solely responsible for declaring that the Medical Devices mentioned in this statement are of Low-Risk Class (Class I) and comply with the requirements of the European Regulation 745/2017 and where appropriate, the standards and legislation referred to.

17. Terms of Warranty:

The product is guaranteed for a period of two (2) years warranty from the purchasing date.

The warranty apply on factory defects & DOES NOT cover damages resulting from misuse, poor maintenance, modification, overuse, or non-compliance with the product's operating instructions. Also DOES NOT cover the parts that can wear out during use or over time, like the upholstery, wheels, wheel tire, safety belt, footrests or brakes.

Damages or defects caused by: natural disasters, unauthorized maintenance or repair, power supply problems (where applicable) & carriers are not covered by the warranty.

No service or spare part is covered unless it has been inspected by an authorized service department or the manufacturer of the product. Shipping costs for goods of warranty are to be paid by the user- customer or by dealer.

When the user or the provider apply on warranty request, provide a valid proof of purchase and the completed warranty form, with a photo of the product serial number and a photo or video of the problem that has occurred. In the event the user cannot provide the above documents, if we can confirm that the product or part is still within the warranty period, we may be able to honor the user's request for coverage under warranty, but we are not obligated to do so.

18. Warranty Form

`BUYERS DETAILS:	
FULL NAME:	
PHONE:	
PURCHASING DATE:	
INVOICE NUMBER:	
SERIAL NUMBER:	
LOT:	
DEALER DETAILS:	
FULL NAME:	
ADDRESS:	
PHONE:	
	DEALER SIGN AND SIGNATURE

MOBIAK S.A
ΚΑΘΙΑΝΑ ΑΚΡΩΤΗΡΙΟΥ, ΧΑΝΙΑ - ΚΡΗΤΗ
KATHIANA AKROTIRIOU CHANIA - CRETE
T: +30 28210 63222 | WEB: www.mobiakcare.com

CE

MOBIAK®
MEDICAL • GASES • FIREFIGHTING est. 1977

19. Disclaimer:

The content of the product manual is the intellectual property of MOBIAK S.A. (hereinafter referred to as the "Company") and the copyright is reserved by the Company. Without the written consent of the Company, no one may copy, export or translate the manual or any content included in any form into other languages without the written consent of the Company. The Company carefully compiles the product manual with a sense of responsibility for users but does not guarantee the completeness or correctness of its contents. The product manual, as a purely technical document, does not include implied content or content with oblique reference to third parties. The Company assumes no responsibility for any typographical errors that may cause ambiguities. In the event of direct or indirect loss of information or business interruption due to the product manual or the product information contained therein, the Company shall not be liable.

When using all MOBIAK S.A. products, users must strictly follow the product manual and the instructions for use available on the Company's official website (www.mobiakcare.com) for the handling of the equipment and not violate any warning marked in the product manual when using the equipment. In case of injury to users or other persons, or material damage caused by the handling of the equipment without following the instructions for use or by violating the warnings marked in the product manual, the persons concerned shall bear the respective legal responsibilities.

As the Company continues to improve and update its products, it is essential that you regularly log on to its official website to refer to the most up-to-date instructions and conditions of use or contact your distributor for the latest user manual. The Company reserves the right to revise the product manual without notice.

The Company reserves the right of final interpretation.

MOBIAK S.A
ΚΑΘΙΑΝΑ ΑΚΡΩΤΗΡΙΟΥ, ΧΑΝΙΑ - ΚΡΗΤΗ
KATHIANA AKROTIRIOU CHANIA - CRETE
T: +30 28210 63222 | WEB: www.mobiakcare.com

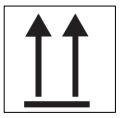
CE

MOBIAK® MEDICAL • GASES • FIREFIGHTING est. 1977

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Description of symbols that have been printed on the label, Master Carton and or User Manual

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